


DHEA			in pg/mL		
Lab Code	Method	Producer	Sample 1	Sample 2	Sample 3
3013	LUM	IBL	474,8	99,1	286,1
3033	ELISA	IBL	518,7	121,0	331,4
3047	GC/LC-MS	Others	371,8	46,0	196,8
3048	ELISA	IBL	449,6	66,3	288,2
3060	ELISA	IBL	411,5	109,3	308,4
3061	ELISA	IBL	464,7	91,3	361,2
13004	ELISA	IBL	423,9	101,2	287,4
15033	GC/LC-MS	Others	429,4	93,7	250,8
20003	ELISA	Others	410,3	77,8	267,6
22236	ELISA	Salimetrics	489,0	74,0	301,0
30001	LUM	IBL	547,6	135,4	426,5
40001	ELISA	IBL	547,0	142,0	367,0
40002	ELISA	IBL	559,9	115,3	337,1
40003	GC/LC-MS	Others	621,1	229,0	641,6
40013	ELISA	DRG	730,0	188,0	460,0
40015	ELISA	Demeditec	668,7	132,3	493,6
40016	ELISA	Demeditec	689,0	123,3	386,5
40017	ELISA	DRG	830,5	159,2	551,0
40024	ELISA	IBL	502,1	108,5	330,5
40029	ELISA	Demeditec	722,9	141,6	424,8
40036	ELISA	IBL	441,7	105,4	288,7
40041	ELISA	Demeditec	669,4	106,0	371,3
40049	ELISA	IBL	465,7	102,7	302,6
40055	ELISA	IBL	460,0	109,0	290,0


Excluded Samples (as probably reported for DHEA-S)

			in pg/mL		
Lab Code	Producer		Sample 1	Sample 2	Sample 3
3051	ELISA	IBL	3300,00	1700,00	1500,00

Sample 1								
Method		LUM	ELISA	ELISA	ELISA	ELISA	ELISA	GC/LC-MS
Supplier	all	IBL	IBL	Salimetrics	Demeditec	DRG	Others	Others
Mean	537,5	511,2	476,8	489,0	687,5	780,3	410,3	474,1
SD	120,4	51,4	48,8		25,4	71,1		130,5
Mean - 2SD	296,7	408,3	379,2		636,7	638,1		213,1
Mean + 2SD	778,3	614,1	574,4		738,3	922,4		735,1
Min.	371,8	474,8	411,5	489,0	668,7	730,0	410,3	371,8
Max.	830,5	547,6	559,9	489,0	722,9	830,5	410,3	621,1
Amount	24	2	11	1	4	2	1	3

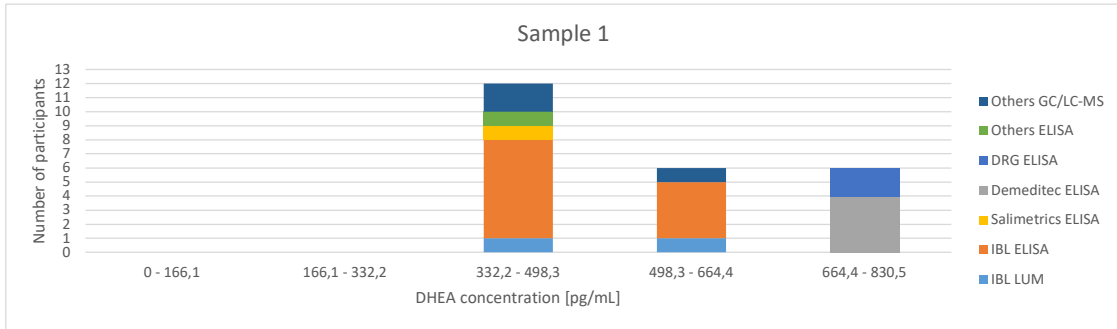

Sample 2								
Method		LUM	ELISA	ELISA	ELISA	ELISA	ELISA	GC/LC-MS
Supplier	all	IBL	IBL	Salimetrics	Demeditec	DRG	Others	Others
Mean	115,7	117,2	106,5	74,0	125,8	173,6	77,8	122,9
SD	38,0	25,7	18,6		15,2	20,4		94,9
Mean - 2SD	39,8	65,8	69,3		95,5	132,9		0,0
Mean + 2SD	191,7	168,7	143,8		156,1	214,3		312,8
Min.	46,0	99,1	66,3	74,0	106,0	159,2	77,8	46,0
Max.	229,0	135,4	142,0	74,0	141,6	188,0	77,8	229,0
Amount	24	2	11	1	4	2	1	3

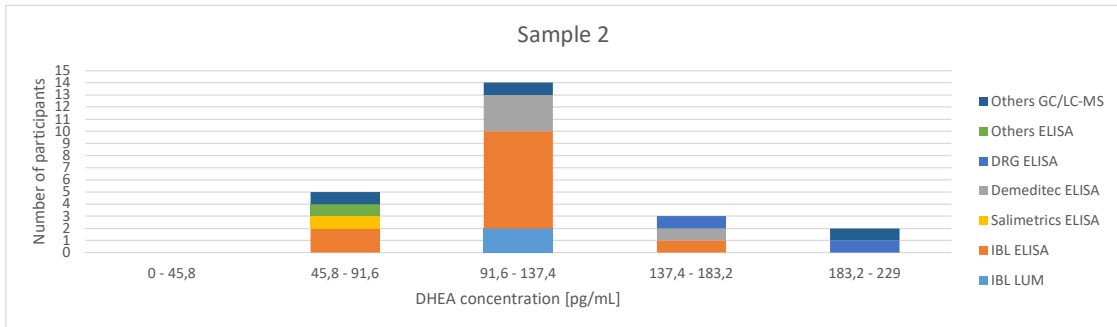

Sample 3								
Method		LUM	ELISA	ELISA	ELISA	ELISA	ELISA	GC/LC-MS
Supplier	all	IBL	IBL	Salimetrics	Demeditec	DRG	Others	Others
Mean	356,3	356,3	317,5	301,0	419,0	505,5	267,6	363,1
SD	99,5	99,3	29,6		54,6	64,3		242,7
Mean - 2SD	157,3	157,8	258,3		309,9	376,8		0,0
Mean + 2SD	555,2	554,8	376,7		528,2	634,2		848,5
Min.	196,8	286,1	287,4	301,0	371,3	460,0	267,6	196,8
Max.	641,6	426,5	367,0	301,0	493,6	551,0	267,6	641,6
Amount	24	2	11	1	4	2	1	3

