



Progesterone			in pg/mL			
Lab Code	Method	Producer	Sample 1	Sample 2	Sample 3	
3013	LUM	IBL	63,69	17,89	219,95	
3033	ELISA	IBL	88,33	50,47	252,37	
3047	GC/LC-MS	Others	61,81	14,71	200,75	
3051	ELISA	IBL	59,40	2,70	220,30	
3060	ELISA	IBL	76,34	19,52	209,62	
3061	ELISA	IBL	74,74	32,55	236,30	
13004	LUM	IBL	95,20	22,20	255,70	
15033	GC/LC-MS	Others	77,20	6,30	266,40	
20002	ELISA	IBL	51,37	17,81	224,69	
20003	ELISA	Others	76,20	36,50	262,90	
22236	ELISA	Salimetrics	94,32	7,26	289,27	
30001	LUM	IBL	81,61	17,67	306,75	
40001	ELISA	IBL	89,00	28,00	272,00	
40002	ELISA	IBL	86,20	37,75	291,15	
40003	GC/LC-MS	Others	56,87	16,36	198,08	
40013	ELISA	DRG	108,82	18,85	390,49	
40014	GC/LC-MS	Others	31,00	0,00	190,00	
40015	ELISA	Demeditec	130,00	17,70	439,00	
40016	ELISA	DRG	128,10	18,80	371,00	
40017	ELISA	DRG	154,80	40,20	491,60	
40024	ELISA	IBL	82,66	31,15	266,30	
40036	ELISA	IBL	71,40	26,70	209,67	
40041	ELISA	IBL	85,10	28,90	242,00	
40049	ELISA	IBL	87,41	32,27	200,28	
40055	ELISA	IBL	99,80	38,40	285,00	
40057	ELISA	IBL	117,10	61,60	333,90	

Sample 1									
Method		LUM	ELISA	ELISA	ELISA	ELISA	ELISA	GC/LC-MS	
Supplier	all	IBL	IBL	Salimetrics	Demeditec	DRG	Others	Others	
Mean	85,71	80,17	82,22	94,32	130,00	130,57	76,20	56,72	
SD	26,16	15,80	16,68			23,09		19,21	
Mean - 2SD	33,40	48,56	48,85			84,39		18,30	
Mean + 2SD	138,02	111,78	115,59			176,75		95,14	
Min.	31,00	63,69	51,37	94,32	130,00	108,82	76,20	31,00	
Max.	154,80	95,20	117,10	94,32	130,00	154,80	76,20	77,20	
Amount	26	3	13	1	1	3	1	4	

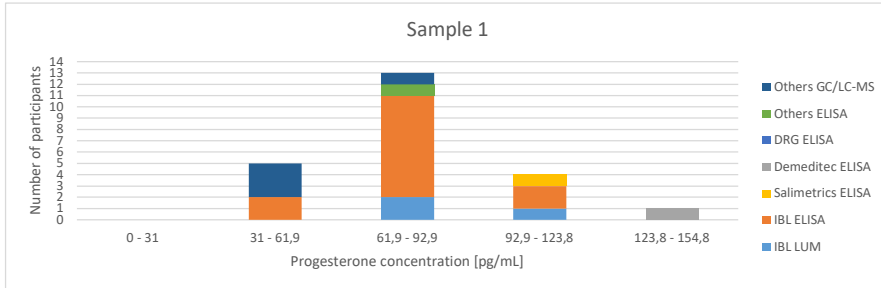

Sample 2									
Method		LUM	ELISA	ELISA	ELISA	ELISA	ELISA	GC/LC-MS	
Supplier	all	IBL	IBL	Salimetrics	Demeditec	DRG	Others	Others	
Mean	24,70	19,25	31,37	7,26	17,70	25,95	36,50	9,34	
SD	14,08	2,56	14,60			12,34		7,63	
Mean - 2SD	0,00	14,14	2,18			1,27		0,00	
Mean + 2SD	52,87	24,36	60,56			50,63		24,60	
Min.	0,00	17,67	2,70	7,26	17,70	18,80	36,50	0,00	
Max.	61,60	22,20	61,60	7,26	17,70	40,20	36,50	16,36	
Amount	26	3	13	1	1	3	1	4	

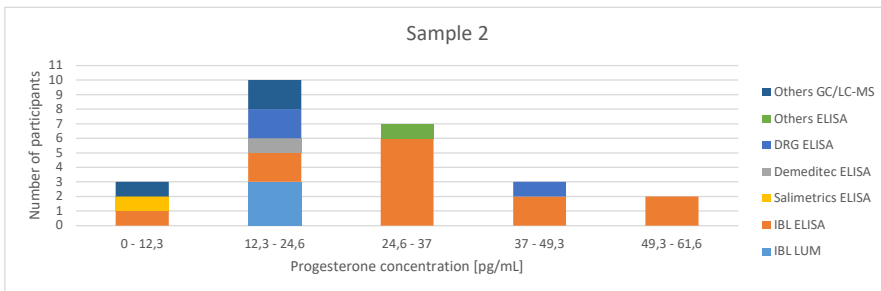

Sample 3									
Method		LUM	ELISA	ELISA	ELISA	ELISA	ELISA	GC/LC-MS	
Supplier	all	IBL	IBL	Salimetrics	Demeditec	DRG	Others	Others	
Mean	274,06	260,80	249,51	289,27	439,00	417,70	262,90	213,81	
SD	75,06	43,62	38,96			64,74		35,36	
Mean - 2SD	123,93	173,55	171,59			288,22		143,09	
Mean + 2SD	424,18	348,05	327,42			547,18		284,52	
Min.	190,00	219,95	200,28	289,27	439,00	371,00	262,90	190,00	
Max.	491,60	306,75	333,90	289,27	439,00	491,60	262,90	266,40	
Amount	26	3	13	1	1	3	1	4	

