



Testosterone			in pg/mL		
Lab Code	Method	Producer	Sample 1	Sample 2	Sample 3
3013	LUM	IBL	194,04	36,02	11,32
3019	GC/LC-MS	Others	185,22	36,66	8,65
3030	GC/LC-MS	Others	198,24	39,31	9,65
3032	GC/LC-MS	Others	127,67	26,22	7,20
3033	ELISA	IBL	213,26	46,11	24,50
3047	GC/LC-MS	Others	204,07	46,17	10,73
3051	ELISA	IBL	235,00	37,70	16,00
3060	ELISA	IBL	130,10	16,54	0,92
3061	ELISA	IBL	233,91	62,23	25,91
13004	LUM	IBL	202,60	33,40	9,40
15033	GC/LC-MS	Others	224,50	43,50	10,90
20002	ELISA	IBL	238,88	39,73	15,12
20003	ELISA	Others	209,60	54,30	30,40
22236	LUM	IBL	266,28	51,59	16,14
30001	LUM	IBL	227,52	38,65	12,16
40001	LUM	IBL	210,20	44,50	17,50
40002	ELISA	IBL	239,34	50,81	21,09
40003	GC/LC-MS	Others	182,02	44,97	13,81
40013	ELISA	DRG	230,00	52,30	18,80
40014	GC/LC-MS	Others	210,00	45,00	7,00
40015	ELISA	Demeditec	201,00	30,40	11,70
40016	ELISA	Demeditec	199,30	55,40	33,40
40017	ELISA	DRG	279,80	41,70	9,80
40024	LUM	IBL	260,30	37,00	13,30
40029	ELISA	Demeditec	186,21	56,11	33,71
40036	ELISA	IBL	210,60	42,70	24,40
40041	ELISA	Demeditec	277,90	63,60	37,80
40049	ELISA	IBL	234,04	44,10	25,02
40054	GC/LC-MS	Others	138,00	30,00	5,00
40055	ELISA	IBL	228,00	59,00	31,20
40057	ELISA	IBL	199,60	47,60	24,40

Sample 1							
Method		LUM	ELISA	ELISA	ELISA	ELISA	GC/LC-MS
Supplier	all	IBL	IBL	Demeditec	DRG	Others	Others
Mean	212,2	226,8	216,3	216,1	254,9	209,6	183,7
SD	36,3	30,4	33,2	41,7	35,2		34,3
Mean - 2SD	139,6	166,1	149,9	132,7	184,5		115,2
Mean + 2SD	284,8	287,6	282,6	299,6	325,3		252,2
Min.	127,7	194,0	130,1	186,2	230,0	209,6	127,7
Max.	279,8	266,3	239,3	277,9	279,8	209,6	224,5
Amount	31	6	10	4	2	1	8

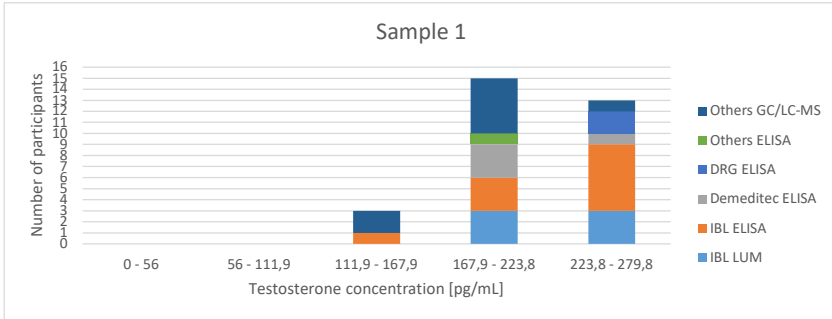

Sample 2							
Method		LUM	ELISA	ELISA	ELISA	ELISA	GC/LC-MS
Supplier	all	IBL	IBL	Demeditec	DRG	Others	Others
Mean	43,65	40,19	44,65	51,38	47,00	54,30	38,98
SD	10,42	6,70	12,61	14,47	7,50		7,49
Mean - 2SD	22,82	26,79	19,43	22,44	32,01		23,99
Mean + 2SD	64,49	53,59	69,87	80,31	61,99		53,97
Min.	16,54	33,40	16,54	30,40	41,70	54,30	26,22
Max.	63,60	51,59	62,23	63,60	52,30	54,30	46,17
Amount	31	6	10	4	2	1	8

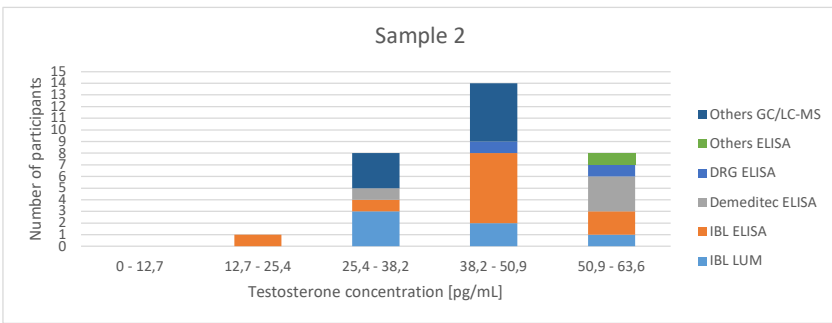

Sample 3							
Method		LUM	ELISA	ELISA	ELISA	ELISA	GC/LC-MS
Supplier	all	IBL	IBL	Demeditec	DRG	Others	Others
Mean	17,32	13,30	20,86	29,15	14,30	30,40	9,12
SD	9,35	3,04	8,45	11,81	6,36		2,76
Mean - 2SD	0,00	7,23	3,96	5,54	1,57		3,59
Mean + 2SD	36,02	19,38	37,75	52,77	27,03		14,64
Min.	0,92	9,40	0,92	11,70	9,80	30,40	5,00
Max.	37,80	17,50	31,20	37,80	18,80	30,40	13,81
Amount	31	6	10	4	2	1	8

