


DHEA			in pg/mL		
Lab Code	Method	Producer	Sample 1	Sample 2	Sample 3
3013	LUM	IBL	452,7	117,1	322,4
3033	ELISA	IBL	533,1	121,0	289,0
3048	ELISA	IBL	446,7	121,0	371,8
3060	ELISA	IBL	428,1	94,5	262,0
3061	ELISA	IBL	373,3	80,1	179,7
13004	ELISA	IBL	504,4	92,7	250,8
15033	GC/LC-MS	Others	357,4	67,2	272,7
22236	ELISA	Salimetrics	562,0	54,0	321,0
30001	LUM	IBL	417,9	121,0	334,3
40001	ELISA	IBL	597,0	137,0	368,0
40002	ELISA	IBL	539,8	117,0	338,8
40013	ELISA	DRG	768,0	166,0	471,0
40015	ELISA	Demeditec	1180,0	144,3	480,5
40016	ELISA	Demeditec	229,7	38,9	128,8
40017	ELISA	DRG	767,0	106,2	466,5
40024	ELISA	IBL	495,1	112,4	327,8
40027	ELISA	IBL	527,1	116,8	350,9
40029	ELISA	Demeditec	797,5	151,6	458,7
40036	ELISA	IBL	475,2	111,9	312,5
40041	ELISA	Demeditec	641,4	117,3	435,2
40049	ELISA	IBL	423,6	104,7	258,6
40055	ELISA	IBL	546,0	103,0	302,0
60006	GC/LC-MS	Others	464,0	85,6	229,1


Excluded Samples (as probably reported for DHEA-S)

			in pg/mL		
Lab Code	Method	Producer	Sample 1	Sample 2	Sample 3
3036	ELISA	DRG	2391,9	1959,7	2161,4
13009	ELISA	IBL	2300,0	1500,0	2100,0

Sample 1							
Method		LUM	ELISA	ELISA	ELISA	ELISA	GC/LC-MS
Supplier	all	IBL	IBL	Salimetrics	Demeditec	DRG	Others
Mean	544,6	435,3	490,8	562,0	712,2	767,5	410,7
SD	189,0	24,6	63,6		393,3	0,7	75,4
Mean - 2SD	166,7	386,1	363,6		0,0	766,1	260,0
Mean + 2SD	922,6	484,5	617,9		1498,7	768,9	561,4
Min.	229,7	417,9	373,3	562,0	229,7	767,0	357,4
Max.	1180,0	452,7	597,0	562,0	1180,0	768,0	464,0
Amount	23	2	12	1	4	2	2

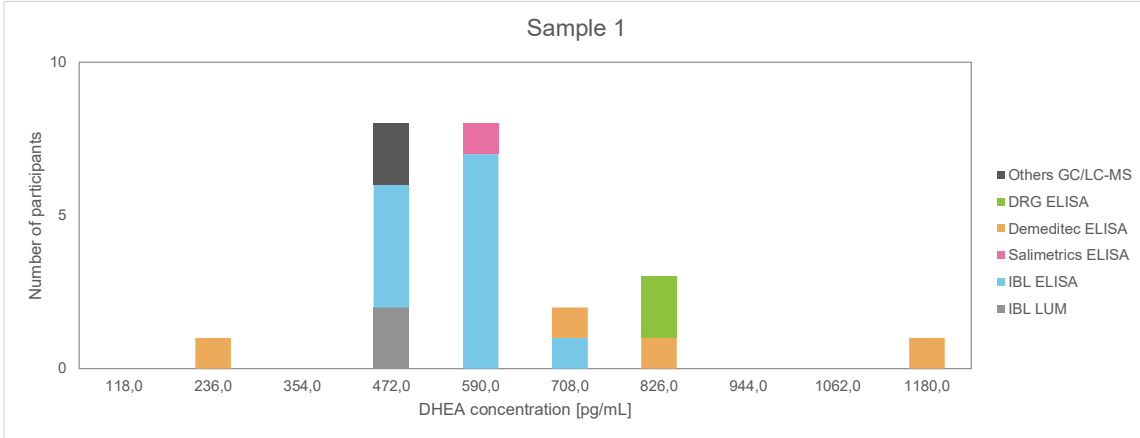

Sample 2							
Method		LUM	ELISA	ELISA	ELISA	ELISA	GC/LC-MS
Supplier	all	IBL	IBL	Salimetrics	Demeditec	DRG	Others
Mean	107,9	119,1	109,3	54,0	113,0	136,1	76,4
SD	29,0	2,8	15,3		51,6	42,3	13,0
Mean - 2SD	49,8	113,5	78,7		9,9	51,5	50,4
Mean + 2SD	165,9	124,7	140,0		216,2	220,7	102,4
Min.	38,9	117,1	80,1	54,0	38,9	106,2	67,2
Max.	166,0	121,0	137,0	54,0	151,6	166,0	85,6
Amount	23	2	12	1	4	2	2

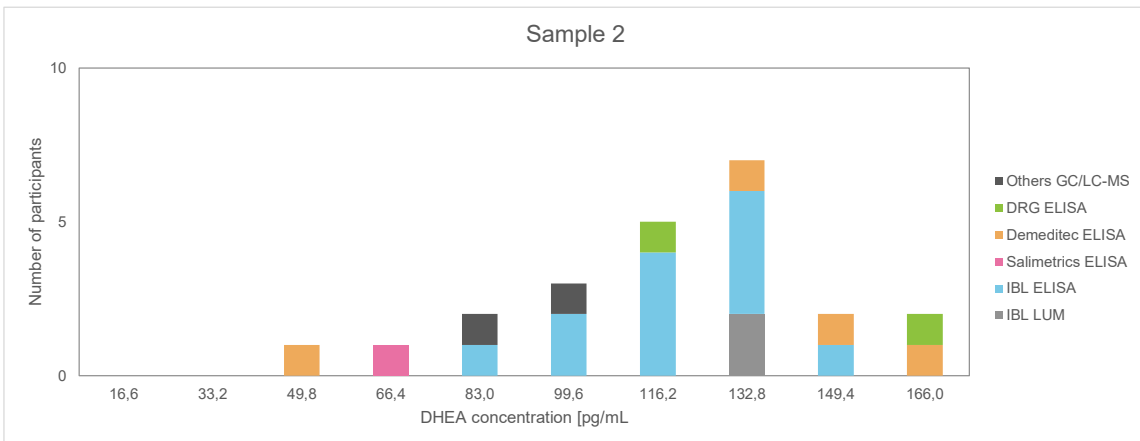

Sample 3							
Method		LUM	ELISA	ELISA	ELISA	ELISA	GC/LC-MS
Supplier	all	IBL	IBL	Salimetrics	Demeditec	DRG	Others
Mean	327,5	328,3	301,0	321,0	375,8	468,8	250,9
SD	90,4	8,4	56,4		165,7	3,2	30,8
Mean - 2SD	146,8	311,4	188,2		44,4	462,4	189,3
Mean + 2SD	508,2	345,2	413,8		707,2	475,1	312,6
Min.	128,8	322,4	179,7	321,0	128,8	466,5	229,1
Max.	480,5	334,3	371,8	321,0	480,5	471,0	272,7
Amount	23	2	12	1	4	2	2

