



Progesterone			in pg/mL		
Lab Code	Method	Producer	Sample 1	Sample 2	Sample 3
3013	LUM	IBL	82,1	20,1	270,2
3033	ELISA	IBL	116,7	44,2	340,7
3036	ELISA	DRG	141,9	16,7	428,9
3051	ELISA	IBL	133,0	35,8	266,4
3060	ELISA	IBL	83,1	33,6	261,8
3061	ELISA	IBL	59,9	24,4	163,7
13004	LUM	IBL	75,0	20,1	264,6
13009	ELISA	IBL	148,7	45,0	394,1
15033	GC/LC-MS	Others	82,6	6,6	296,7
22236	ELISA	Others	86,8	6,3	358,7
30001	LUM	IBL	43,2	13,7	294,0
40001	ELISA	IBL	80,0	36,0	280,0
40002	ELISA	IBL	77,0	29,8	255,7
40013	ELISA	DRG	151,0	24,1	440,0
40015	ELISA	Others	118,0	12,4	451,0
40016	ELISA	DRG	42,6	6,3	146,8
40017	ELISA	DRG	133,0	16,6	481,9
40024	ELISA	IBL	86,2	32,3	256,3
40027	ELISA	IBL	101,2	33,2	293,1
40036	ELISA	IBL	75,7	33,3	251,7
40041	ELISA	IBL	69,9	25,9	284,2
40049	ELISA	IBL	81,8	30,5	303,3
40055	ELISA	IBL	90,0	28,9	235,0
2002	ELISA	IBL	94,7	31,1	289,4

Sample 1						
Method		LUM	ELISA	ELISA	ELISA	GC/LC-MS
Supplier	all	IBL	IBL	DRG	Others	Others
Mean	93,9	66,8	92,7	117,1	102,4	82,6
SD	30,0	20,7	24,8	50,2	22,1	
Mean - 2SD	34,0	25,4	43,2	16,7	58,2	
Mean + 2SD	153,8	108,2	142,2	217,6	146,6	
Min.	42,6	43,2	59,9	42,6	86,8	82,6
Max.	151,0	82,1	148,7	151,0	118,0	82,6
Amount	24	3	14	4	2	1

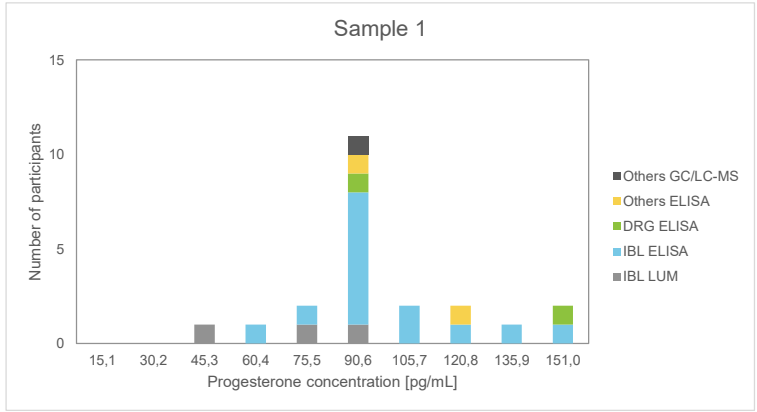

Sample 2						
Method		LUM	ELISA	ELISA	ELISA	GC/LC-MS
Supplier	all	IBL	IBL	DRG	Others	Others
Mean	25,3	18,0	33,1	15,9	9,4	6,6
SD	11,0	3,7	5,9	7,3	4,3	
Mean - 2SD	3,4	10,6	21,4	1,3	0,7	
Mean + 2SD	47,2	25,3	44,9	30,5	18,0	
Min.	6,3	13,7	24,4	6,3	6,3	6,6
Max.	45,0	20,1	45,0	24,1	12,4	6,6
Amount	24	3	14	4	2	1

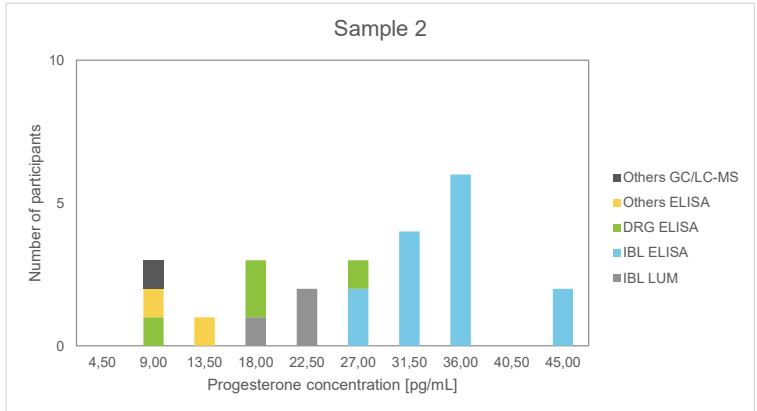

Sample 3						
Method		LUM	ELISA	ELISA	ELISA	GC/LC-MS
Supplier	all	IBL	IBL	DRG	Others	Others
Mean	304,5	276,3	276,8	374,4	404,8	296,7
SD	82,6	15,6	52,2	153,4	65,3	
Mean - 2SD	139,4	245,0	172,4	67,6	274,3	
Mean + 2SD	469,6	307,5	381,2	681,2	535,4	
Min.	146,8	264,6	163,7	146,8	358,7	296,7
Max.	481,9	294,0	394,1	481,9	451,0	296,7
Amount	24	3	14	4	2	1

