

## Quality Assessment Scheme for Salivary Steroids

### Estradiol Sample Number 23-21


Your result: The results from all participants can be found in the outstanding table.

Method		LUM	ELISA		MS	
Supplier	all	IBL	IBL	Demeditec	DRG	In-House
Lab. Code	pg/mL	pg/mL	pg/mL	pg/mL	pg/mL	pg/mL
3013	12,70	12,70				
3036	0,865				0,865	
3060	7,788		7,788			
3061	3,700		3,700			
6005	4,618		4,618			
13004	11,50	11,50				
15033	7,600					7,600
20002	7,140		7,140			
22236	11,343	11,34				
30001	11,226	11,23				
40001	7,400		7,400			
40002	8,560		8,560			
40003-1	6,200					6,200
40013	3,900				3,900	
40015	2,920			2,920		
40016	5,100				5,100	
40024	9,200	9,200				
40027	5,280		5,280			
40036	7,850		7,850			
40041	3,050			3,050		
40049	8,763		8,763			
40055	7,520		7,520			
50002	8,400		8,400			

### Statistical Results

Method		LUM	ELISA		MS	
Supplier	all	IBL	IBL	Demeditec	DRG	In-House
Mean	7,071	11,194	7,002	2,985	3,288	6,900
SD	3,068	1,261	1,698	0,092	2,183	0,990
Mean - 2SD	0,934	8,672	3,606	2,801	0,000	4,920
Mean + 2SD	13,207	13,716	10,397	3,169	7,654	8,880
Amount	23	5	11	2	3	2

\*The graphic is for information purposes only and not validated. Please read your results from the table.

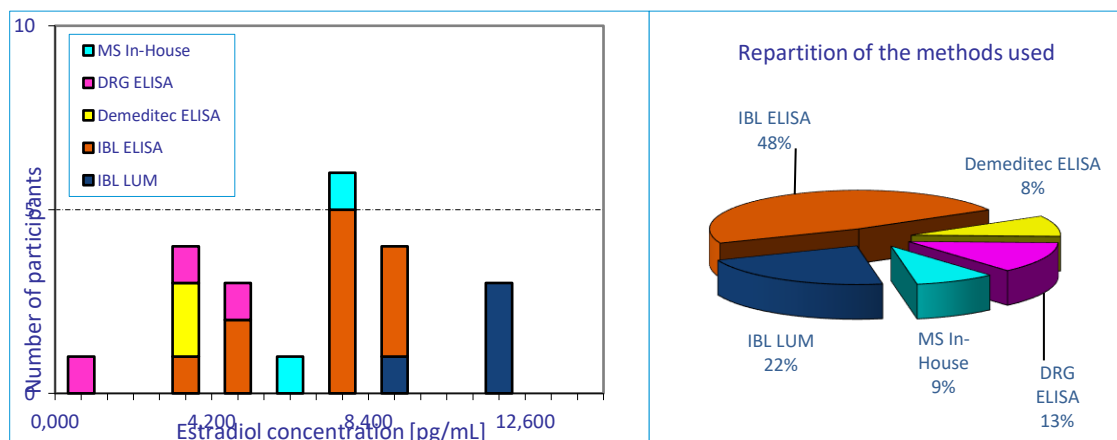

## Quality Assessment Scheme for Salivary Steroids

### Estradiol Sample Number 23-22


Your result: The results from all participants can be found in the outstanding table.

Method		LUM	ELISA		MS
Supplier	all	IBL	IBL	Demeditec	In-House
Lab. Code	pg/mL	pg/mL	pg/mL	pg/mL	pg/mL
3013	2,280	2,280			
3036	0,836				0,836
3060	1,421		1,421		
3061	2,250		2,250		
6005	1,873		1,873		
13004	3,900	3,900			
15033	0,300				0,300
20002	2,640		2,640		
22236	2,095	2,095			
30001	2,330	2,330			
40001	2,500		2,500		
40002	3,140		3,140		
40003-1	1,340				1,340
40013	2,870				2,870
40015	1,990			1,990	
40016	3,400				3,400
40024	0,800	0,800			
40027	2,150		2,150		
40036	2,290		2,290		
40041	1,600			1,600	
40049	2,704		2,704		
40055	2,330		2,330		
50002	2,100		2,100		

### Statistical Results

Method		LUM	ELISA		MS
Supplier	all	IBL	IBL	Demeditec	In-House
Mean	2,136	2,281	2,309	1,795	0,820
SD	0,839	1,101	0,452	0,276	0,735
Mean - 2SD	0,459	0,078	1,405	1,243	0,000
Mean + 2SD	3,814	4,484	3,213	2,347	2,291
Amount	23	5	11	2	3

\*The graphic is for information purposes only and not validated. Please read your results from the table.

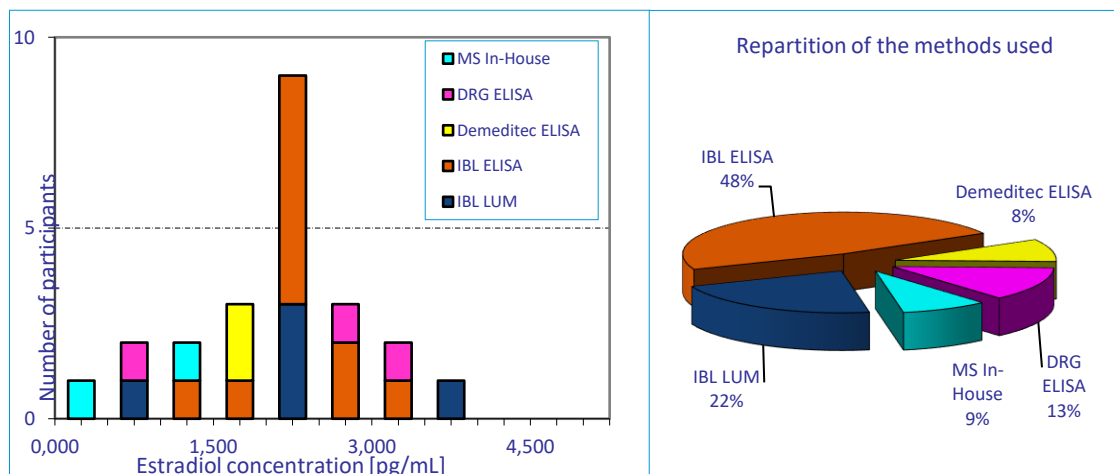

## Quality Assessment Scheme for Salivary Steroids

### Estradiol Sample Number 23-23

Your result: The results from all participants can be found in the outstanding table.

Method		LUM	ELISA		MS
Supplier	all	IBL	IBL	Demeditec	In-House
Lab. Code	pg/mL	pg/mL	pg/mL	pg/mL	pg/mL
3013	5,970	5,970			
3036	2,709				2,709
3060	3,543		3,543		
3061	3,930		3,930		
6005	2,880		2,880		
13004	6,300	6,300			
15033	2,500				2,500
20002	4,110		4,110		
22236	5,755	5,755			
30001	4,932	4,932			
40001	4,200		4,200		
40002	4,530		4,530		
40003-1	3,100				3,100
40013	3,450				3,450
40015	2,530			2,530	
40016	4,700				4,700
40024	4,000	4,000			
40027	3,390		3,390		
40036	4,310		4,310		
40041	2,250			2,250	
40049	5,528		5,528		
40055	4,450		4,450		
50002	4,200		4,200		

### Statistical Results

Method		LUM	ELISA		MS
Supplier	all	IBL	IBL	Demeditec	In-House
Mean	4,055	5,391	4,097	2,390	3,620
SD	1,141	0,927	0,689	0,198	1,006
Mean - 2SD	1,773	3,537	2,720	1,994	1,607
Mean + 2SD	6,337	7,246	5,475	2,786	5,632
Amount	23	5	11	2	3

\*The graphic is for information purposes only and not validated. Please read your results from the table.

