

## Quality Assessment Scheme for Salivary Steroids

### Progesterone Sample Number 24-11


Your result: The results from all participants can be found in the outstanding table.

Method		LUM	ELISA		MS
Supplier	all	IBL	IBL	DRG	In-House
Lab. Code	pg/mL	pg/mL	pg/mL	pg/mL	pg/mL
3013	75,20	75,20			
3047	49,39				49,39
3061	52,03		52,03		
13004	73,10	73,10			
15033	78,80				78,80
30001	76,66	76,66			
40001	90,00		90,00		
40002	97,85		97,85		
40003-1	51,86				51,86
40014	66,00				66,00
40016	136,00			136,00	
40024	86,00		86,00		
40027	82,61		82,61		
40036	89,35		89,35		
40041	73,20		73,20		
40055	85,80		85,80		

### Statistical Results

Method		LUM	ELISA		MS
Supplier	all	IBL	IBL	DRG	In-House
Mean	78,99	74,99	82,11	136,00	61,51
SD	20,88	1,79	14,02		13,65
Mean - 2SD	37,23	71,41	54,07		34,21
Mean + 2SD	120,75	78,57	110,14		88,82
Amount	16	3	8	1	4

\*The graphic is for information purposes only and not validated. Please read your results from the table.

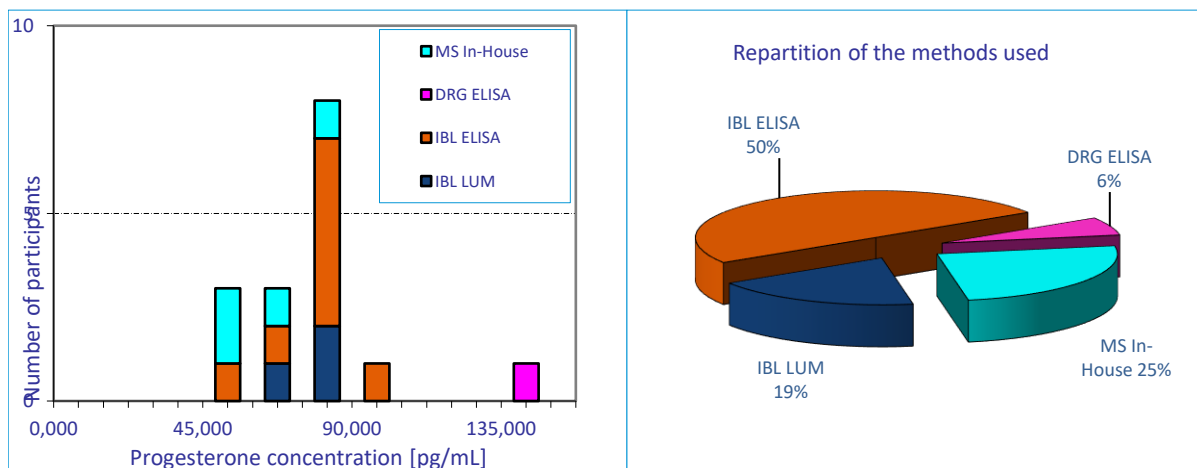

## Quality Assessment Scheme for Salivary Steroids

### Progesterone Sample Number 24-12


Your result: The results from all participants can be found in the outstanding table.

Method		LUM	ELISA		MS
Supplier	all	IBL	IBL	DRG	In-House
Lab. Code	pg/mL	pg/mL	pg/mL	pg/mL	pg/mL
3013	19,10	19,10			
3047	4,03				4,03
3061	29,10		29,10		
13004	19,60	19,60			
15033	6,10				6,10
30001	20,09	20,09			
40001	42,00		42,00		
40002	42,40		42,40		
40003-1	11,98				11,98
40014	5,00				5,00
40016	19,30			19,30	
40024	33,00		33,00		
40027	36,31		36,31		
40036	29,71		29,71		
40041	30,04		30,04		
40055	37,00		37,00		

### Statistical Results

Method		LUM	ELISA		MS
Supplier	all	IBL	IBL	DRG	In-House
Mean	24,05	19,60	34,95	19,30	6,78
SD	12,82	0,50	5,36		3,57
Mean - 2SD	0,000	18,61	24,23		0,000
Mean + 2SD	49,68	20,59	45,66		13,92
Amount	16	3	8	1	4

\*The graphic is for information purposes only and not validated. Please read your results from the table.

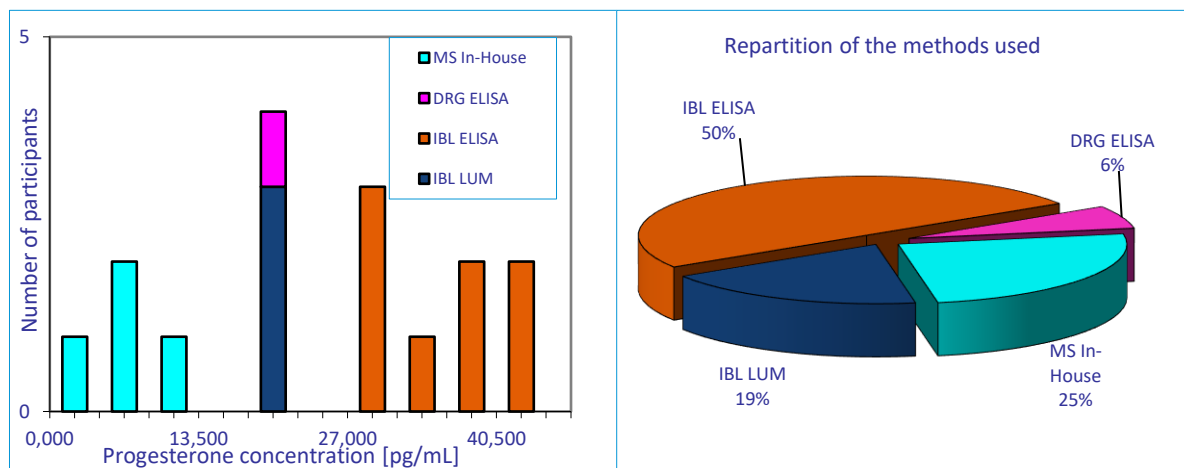

## Quality Assessment Scheme for Salivary Steroids

### Progesterone Sample Number 24-13

Your result: The results from all participants can be found in the outstanding table.

Method		LUM	ELISA		MS
Supplier	all	IBL	IBL	DRG	In-House
Lab. Code	pg/mL	pg/mL	pg/mL	pg/mL	pg/mL
3013	256,00	256,00			
3047	197,72				197,72
3061	219,62		219,62		
13004	243,20	243,20			
15033	296,70				296,70
30001	273,72	273,72			
40001	302,00		302,00		
40002	247,60		247,60		
40003-1	206,08				206,08
40014	265,00				265,00
40016	432,20			432,20	
40024	299,00		299,00		
40027	305,39		305,39		
40036	227,93		227,93		
40041	273,40		273,40		
40055	248,00		248,00		

### Statistical Results

Method		LUM	ELISA		MS
Supplier	all	IBL	IBL	DRG	In-House
Mean	268,35	257,64	265,37	432,20	241,38
SD	55,18	15,33	34,33		47,51
Mean - 2SD	158,00	226,99	196,70		146,36
Mean + 2SD	378,70	288,29	334,03		336,39
Amount	16	3	8	1	4

\*The graphic is for information purposes only and not validated. Please read your results from the table.

