

Quality Assessment Scheme for Salivary Steroids

Testosterone Sample Number 24-11


Your result: The results from all participants can be found in the outstanding table.

Method		LUM	ELISA			MS
Supplier	all	IBL	IBL	Demeditec	DRG	In-House
Lab. Code	pg/mL	pg/mL	pg/mL	pg/mL	pg/mL	pg/mL
3013	200,60	200,60				
3019-1	201,76					201,76
3036	247,26	247,26				
40001	261,70	261,70				
40002	291,63		291,63			
40003-1	237,42					237,42
40014	244,00					244,00
40016	223,50			223,50		
40017	245,40				245,40	
40024	257,00	257,00				
40027	270,27		270,27			
40029	234,27			234,27		
40036	255,03		255,03			
40049	280,80		280,80			
40054	232,00					232,00
40055	250,00		250,00			

Statistical Results

Method		LUM	ELISA			MS
Supplier	all	IBL	IBL	Demeditec	DRG	In-House
Mean	245,79	241,64	269,55	228,89	245,40	228,80
SD	24,88	28,01	17,38	7,62		18,68
Mean - 2SD	196,03	185,61	234,79	213,65		191,44
Mean + 2SD	295,55	297,67	304,30	244,12		266,15
Amount	16	4	5	2	1	4

*The graphic is for information purposes only and not validated. Please read your results from the table.

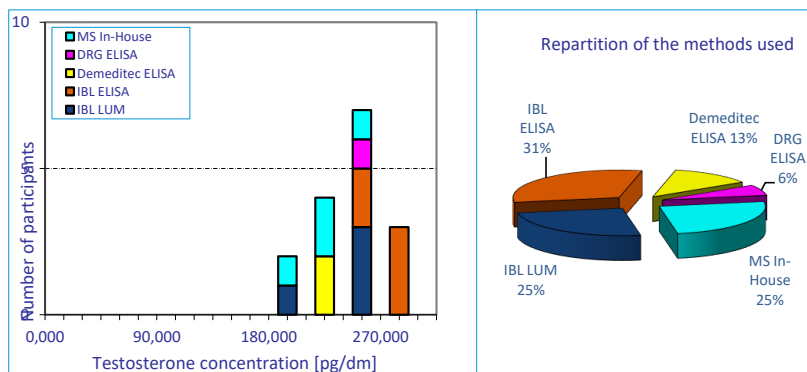

Quality Assessment Scheme for Salivary Steroids

Testosterone Sample Number 24-12


Your result: The results from all participants can be found in the outstanding table.

Method	LUM		ELISA			MS
Supplier	all	IBL	IBL	Demeditec	DRG	In-House
Lab. Code	pg/mL	pg/mL	pg/mL	pg/mL	pg/mL	pg/mL
3013	41,10	41,10				
3019-1	26,95					26,95
3036	44,38	44,38				
40001	50,30	50,30				
40002	57,46		57,46			
40003-1	58,80					58,80
40014	55,00					55,00
40016	59,80			59,80		
40017	52,90				52,90	
40024	39,00	39,00				
40027	55,93		55,93			
40029	63,81			63,81		
40036	42,98		42,98			
40049	55,80		55,80			
40054	47,00					47,00
40055	43,80		43,80			

Statistical Results

Method	LUM		ELISA			MS
Supplier	all	IBL	IBL	Demeditec	DRG	In-House
Mean	49,69	43,70	51,19	61,81	52,90	46,94
SD	9,59	4,93	7,16	2,84		14,21
Mean - 2SD	30,50	33,84	36,87	56,13		18,52
Mean + 2SD	68,88	53,55	65,51	67,48		75,35
Amount	16	4	5	2	1	4

*The graphic is for information purposes only and not validated. Please read your results from the table.

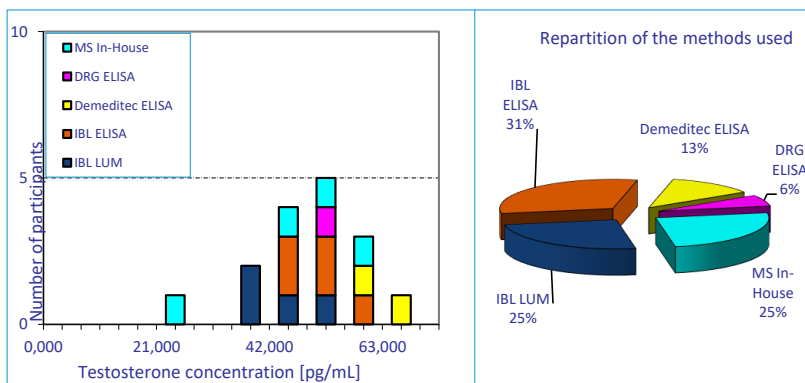

Quality Assessment Scheme for Salivary Steroids

Testosterone Sample Number 24-13

Your result: The results from all participants can be found in the outstanding table.

Method	LUM		ELISA			MS
Supplier	all	IBL	IBL	Demeditec	DRG	In-House
Lab. Code	pg/mL	pg/mL	pg/mL	pg/mL	pg/mL	pg/mL
3013	15,70	15,70				
3019-1	8,65					8,65
3036	19,31	19,31				
40001	12,00	12,00				
40002	24,22		24,22			
40003-1	22,25					22,25
40014	24,00					24,00
40016	30,70			30,70		
40017	25,00				25,00	
40024	14,40	14,40				
40027	22,65		22,65			
40029	37,33			37,33		
40036	21,95		21,95			
40049	25,80		25,80			
40054	26,00					26,00
40055	19,60		19,60			

Statistical Results

Method	LUM		ELISA			MS
Supplier	all	IBL	IBL	Demeditec	DRG	In-House
Mean	21,85	15,35	22,84	34,02	25,00	20,22
SD	7,06	3,05	2,34	4,69		7,87
Mean - 2SD	7,72	9,25	18,15	24,64		4,48
Mean + 2SD	35,97	21,45	27,53	43,39		35,96
Amount	16	4	5	2	1	4

*The graphic is for information purposes only and not validated. Please read your results from the table.

