

Quality Control Certificate



Aldosterone RIA (CT)


REF MG13051	LOT 251705AK	2026-06-27 YYYY-MM-DD
--------------------	---------------------	---------------------------------

Final Release Results:

		Measured value	
Total Activity	TA	41459	cpm
Maximum Binding/TA	B0/TA	62.72	%

	Concentration	Measured value	Range	B/B0
	pg/mL	cpm	cpm	
CAL 0	0.0	26004	-	100%
CAL 1	27.3	21860	-	84%
CAL 2	49.5	18827	-	72%
CAL 3	103.0	15609	-	60%
CAL 4	316.0	10218	-	39%
CAL 5	655.0	7361	-	28%
CAL 6	1512.0	4920	-	19%

	Measured value	Acceptable range		
		lower limit	target value	upper limit
	pg/mL	pg/mL	pg/mL	pg/mL
CONTROL 1	242.53	183.00	261.00	339.00

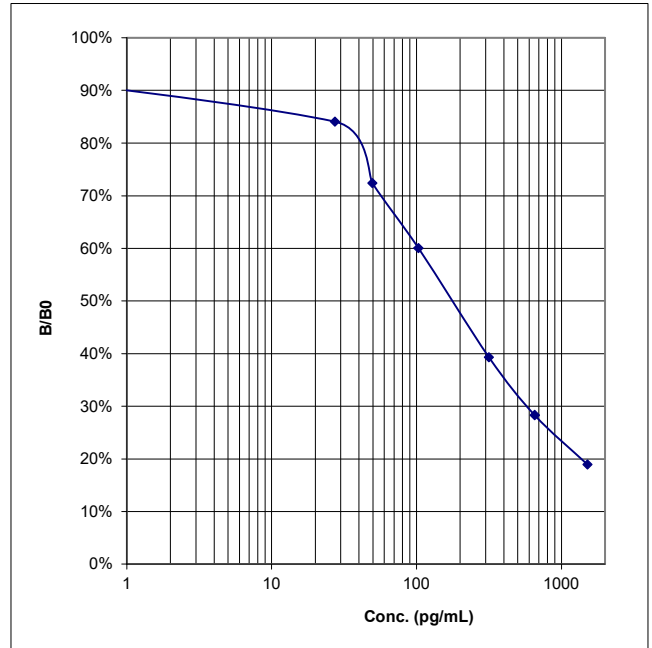

This product has been tested successfully by the Quality Control Department and was released for sale according to the existing specifications.

Verified by: Dr. Spangenberg

Validation date: 2025-04-28

Accepted

This document has been produced electronically and is valid without a signature.



IBL International GmbH
Flughafenstr. 52A, 22335 Hamburg, Germany

Tel.: + 49 (0) 40 532891 -0 Fax: -11
E-MAIL: ibl@tecan.com
WEB: www.tecan.com/ibl