

Quality Control Certificate



17-OH-Progesterone ELISA

REF RE52071

LOT EOH248


2028-07-31
YYYY-MM-DD

Final Release Results:

	Concentration ng/mL	manual			e.g. automated DSX
		Measured value OD	Range* OD	OD/ODmax	Measured value OD
CAL A	0	2.092	> 1.000	100%	2.200
CAL B	0.15	1.524	-	73%	1.465
CAL C	0.5	1.017	-	49%	1.013
CAL D	1.5	0.519	-	25%	0.483
CAL E	3.0	0.344	-	16%	0.325
CAL F	7.5	0.186	-	9%	0.162
CAL G	20	0.092	< 0.350	4%	0.090

*During automated procedure the optical densities of calibrators and controls may rise due to temperature and specific device conditions. In case of none conformance with the OD ranges stated on the QC Certificate please contact IBL-International or use the validation criteria only!

Validation criteria				
	Concentration found ng/mL	Acceptable range		
		lower limit ng/mL	target ng/mL	upper limit ng/mL
CONTROL 1	0.81	0.47	0.86	1.25
CONTROL 2	3.3	2.3	3.5	4.8

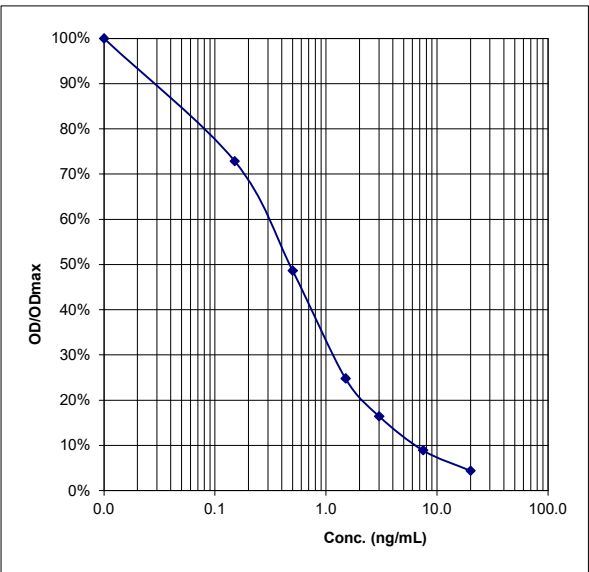

Kit Composition:

REF		LOT	
KEOH651	CONTROL1	K808	2028-07
KEOH661	CONTROL2	K809	2028-07
KEOH511	MTP	K805	2030-03
KEOH541	ENZCONJ	K806	2028-07
KEOH611	CAL A	K807	2028-07
KEOH611	CAL B-G	K807	2028-07
KEZZ751	CONC WASHBUF	K810	2028-07
KEZZ991	TMB SUBS	K610-5	2029-08

REF		LOT	
KEZZ921	TMB STOP	K811	2029-05

Valid edition of the instructions for use: 2024-02

This product has been tested successfully by the Quality Control Department and was released for sale according to the existing specifications.

Verified by: T. Borchert

Validation date: 2026-09-02

Accepted

This document has been produced electronically and is valid without a signature.



IBL International GmbH
Flughafenstr. 52A, 22335 Hamburg, Germany

Tel.: + 49 (0) 40 532891 -0 Fax: -11
E-MAIL: ibl@tecan.com
WEB: www.tecan.com/ibl