

Quality Control Certificate




Amyloid-beta (1-42) CSF ELISA

REF RE59661	LOT EAX129	2026-12-31 YYYY-MM-DD
--------------------	-------------------	---------------------------------

Final Release Results:

	Concentration pg/mL	Measured value OD	Range OD	OD/ODmax
CAL A	0	0.049	< 0.080	2%
CAL B	9.80	0.175		7%
CAL C	19.5	0.306		13%
CAL D	39.0	0.649		27%
CAL E	78.0	1.168		48%
CAL F	156	2.414	-	100%

	Concentration found pg/mL	Acceptable range		
		lower limit pg/mL	target pg/mL	upper limit pg/mL
CONTROL 1	53.4	38.7	55.2	71.8
CONTROL 2	13.7	10.1	14.4	18.7

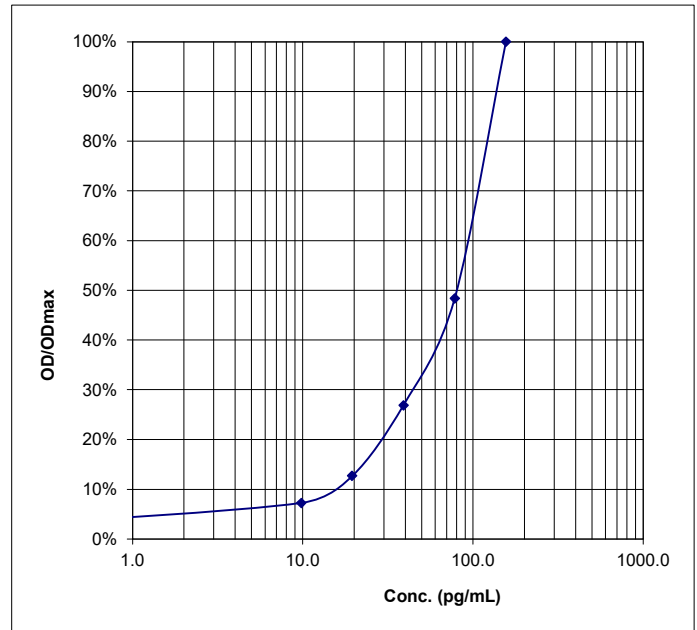

Kit Composition:

REF		LOT		REF		LOT	
KEAX511	MTP	L978	2027-03	KEAX541	CONC ENZCONJ	L979	2027-03
KEAX611	CAL A	L980	2027-03	KEAX571	ASSAYBUF	L874	2026-12
KEAX621	CAL B-F	L981	2027-03	KEZZ451	SAMPLEDIL	L872	2026-12
KEAX651	CONTROL 1	L982	2027-03	KEZZ441	CONC WASHBUF	L938	2027-02
KEAX661	CONTROL 2	L983	2027-03	KEZZ991	TMB SUBS	L541-6	2027-10
				KEZZ921	TMB STOP	L861	2028-01

Valid edition of the instructions for use: 2019-05

This product has been tested successfully by the Quality Control Department and was released for sale according to the existing specifications.

Verified by: Dr. Spangenberg

Validation date: 2025-04-17

Accepted

This document has been produced electronically and is valid without a signature.



IBL International GmbH
Flughafenstr. 52A, 22335 Hamburg, Germany

Tel.: + 49 (0) 40 532891 -0 Fax: -11
E-MAIL: ibl@tecan.com
WEB: www.tecan.com/ibl