

Quality Control Certificate



Progesterone Luminescence Immunoassay

REF RE62021 & RE62029


LOT LPR160

2027-10-31
YYYY-MM-DD

Final Release Results:

	Concentration pg/mL	Measured value RLU	Range RLU	RLU/RLUmax
CAL A	0	1426168	-	100%
CAL B	10	1169938	-	82%
CAL C	25	1031158	-	72%
CAL D	50	818408	-	57%
CAL E	100	638625	-	45%
CAL F	300	302298	-	21%
CAL G	1000	133048	-	9%

	Concentration found pg/mL	Acceptable range		
		lower limit pg/mL	target pg/mL	upper limit pg/mL
CONTROL 1	46.8	31.8	48.9	66.0
CONTROL 2	170	141	177	212

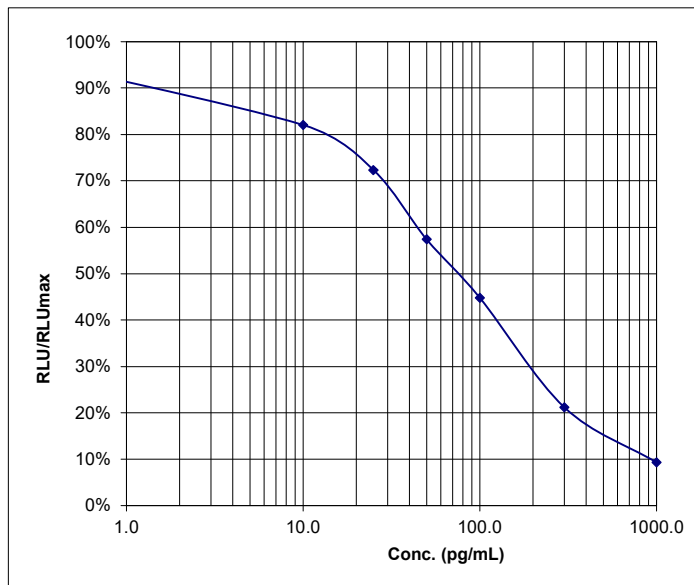

Kit Composition:

REF		LOT	
KLPR651/9	CONTROL 1	K662	2027-10
KLPR661/9	CONTROL 2	K663	2027-10
KLZZ511	MTP	K603	2028-08
KLPR541/9	CONC ENZCONJ	K759	2027-11
KLPR611/9	CAL A	K660	2027-10
KLPR621/9	CAL B-G	K661	2027-10
KLPR521/9	ANTISERUM	K659	2027-10

REF		LOT	
KLZZ681/9	ASSAYBUF	K599	2028-02
KLZZ601/9	LUMINREAG AP	K744	2028-02
KLZZ879	CONC WASHBUF	K586	2028-01

Valid edition of the instructions for use: V2018-07

This product has been tested successfully by the Quality Control Department and was released for sale according to the existing specifications.

Verified by: S. Geweke

Validation date: 2026-07-08

Accepted

This document has been produced electronically and is valid without a signature.



IBL International GmbH
Flughafenstr. 52A, 22335 Hamburg, Germany

Tel.: + 49 (0) 40 532891 -0 Fax: -11
E-MAIL: ibl@tecan.com
WEB: www.tecan.com/ibl