



Progesterone			in pg/mL			
Lab Code	Method	Producer	Sample 1	Sample 2	Sample 3	
3047	GC/LC-MS	Others	86,01	4,09	148,48	
3050	ELISA	IBL	100,04	41,01	243,59	
3065	ELISA	IBL	69,00	34,45	256,30	
3033	ELISA	IBL	75,71	28,39	220,82	
3036	ELISA	DRG	83,85	12,84	327,73	
13004	LUM	IBL	84,70	17,50	262,20	
3061	ELISA	IBL	71,64	36,46	178,52	
13009	ELISA	IBL	106,70	33,04	311,53	
15033	GC/LC-MS	Others	79,00	7,00	255,60	
20002	ELISA	IBL	64,18	25,20	293,42	
20003	ELISA	Others	75,10	16,70	236,90	
22236	ELISA	Salimetrics	105,99	25,87	295,90	
30001	LUM	IBL	47,16	14,79	159,31	
40001	ELISA	IBL	59,00	28,00	228,00	
40002	ELISA	IBL	89,97	36,97	249,03	
40002	ELISA	IBL	89,97	36,97	249,03	
40014	GC/LC-MS	Others	57,00	3,00	138,00	
40015	ELISA	Demeditec	106,70	19,90	291,40	
40016	ELISA	DRG	105,80	16,50	298,00	
40017	ELISA	DRG	158,00	36,10	388,20	
40022	GC/LC-MS	Others	51,10	5,27	179,50	
40024	ELISA	IBL	81,00	30,00	243,00	
40027	ELISA	IBL	99,05	43,24	224,66	
40038	LUM	IBL	67,90	17,95	204,58	
40041	ELISA	IBL	74,40	31,55	250,85	
40049	ELISA	IBL	64,40	31,00	221,60	
40055	ELISA	IBL	66,50	30,60	218,00	

Sample 1								
Method		LUM	ELISA	ELISA	ELISA	ELISA	ELISA	GC/LC-MS
Supplier	all	IBL	IBL	Salimetrics	Demeditec	DRG	Others	Others
Mean	82,22	66,59	79,40	105,99	106,70	115,88	75,10	68,28
SD	22,55	18,80	15,27			38,09		16,85
Mean - 2SD	37,12	28,98	48,86			39,70		34,59
Mean + 2SD	127,31	104,19	109,93			192,06		101,97
Min.	47,16	47,16	59,00	105,99	106,70	83,85	75,10	51,10
Max.	158,00	84,70	106,70	105,99	106,70	158,00	75,10	86,01
Amount	27	3	14	1	1	3	1	4

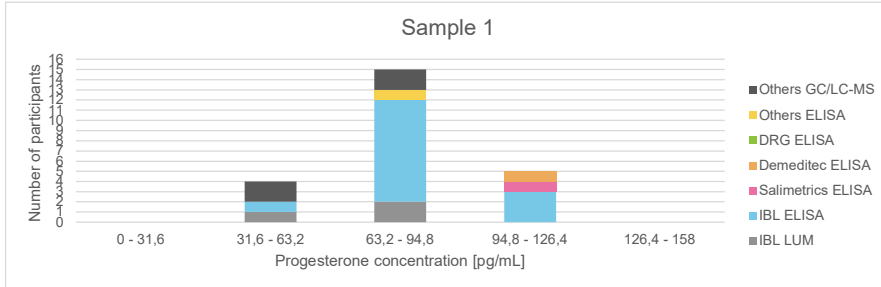

Sample 2								
Method		LUM	ELISA	ELISA	ELISA	ELISA	ELISA	GC/LC-MS
Supplier	all	IBL	IBL	Salimetrics	Demeditec	DRG	Others	Others
Mean	24,61	16,75	33,35	25,87	19,90	21,81	16,70	4,84
SD	11,54	1,71	5,12			12,51		1,71
Mean - 2SD	1,53	13,34	23,11			0,00		1,42
Mean + 2SD	47,69	20,16	43,59			46,83		8,26
Min.	3,00	14,79	25,20	25,87	19,90	12,84	16,70	3,00
Max.	43,24	17,95	43,24	25,87	19,90	36,10	16,70	7,00
Amount	27	3	14	1	1	3	1	4

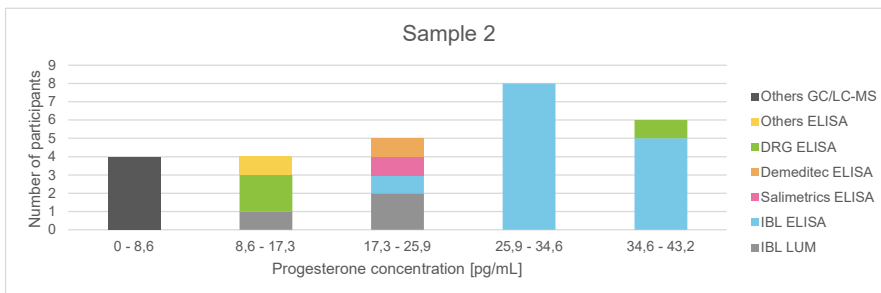

Sample 3								
Method		LUM	ELISA	ELISA	ELISA	ELISA	ELISA	GC/LC-MS
Supplier	all	IBL	IBL	Salimetrics	Demeditec	DRG	Others	Others
Mean	243,49	208,70	242,02	295,90	291,40	337,98	236,90	180,39
SD	55,54	51,57	32,62			45,96		53,14
Mean - 2SD	132,40	105,55	176,79			246,05		74,11
Mean + 2SD	354,58	311,84	307,26			429,91		286,68
Min.	138,00	159,31	178,52	295,90	291,40	298,00	236,90	138,00
Max.	388,20	262,20	311,53	295,90	291,40	388,20	236,90	255,60
Amount	27	3	14	1	1	3	1	4

