



Testosterone

Lab Code	Method	Producer	in pg/mL		
			Sample 1	Sample 2	Sample 3
3019	GC/LC-MS	Others	[862,2]	36,31	8,65
3030	GC/LC-MS	Others	152,19	30,40	9,14
3032	GC/LC-MS	Others	117,00	25,94	4,90
3047	GC/LC-MS	Others	227,01	35,37	9,64
3050	ELISA	IBL	192,43	34,18	13,94
3065	ELISA	IBL	301,50	75,35	33,55
3033	ELISA	IBL	259,37	69,16	31,70
3036	ELISA	DRG	194,81	1,41	1,41
13004	LUM	IBL	178,00	34,00	11,00
3061	ELISA	IBL	209,05	37,66	15,87
13009	ELISA	IBL	306,87	81,18	44,81
15033	GC/LC-MS	Others	213,40	42,70	9,90
20002	ELISA	IBL	245,55	51,45	21,57
20003	ELISA	Others	221,00	61,60	33,70
22236	LUM	IBL	190,49	40,92	7,78
30001	LUM	IBL	214,99	43,80	14,70
40001	LUM	IBL	192,80	39,80	11,20
40002	ELISA	IBL	225,18	39,76	15,60
40002	ELISA	IBL	225,18	39,76	15,60
40014	GC/LC-MS	Others	194,00	30,00	5,00
40015	ELISA	Demeditec	236,40	63,60	40,80
40016	ELISA	Demeditec	223,00	63,90	36,00
40017	ELISA	DRG	293,20	48,30	12,30
40022	GC/LC-MS	Others	160,52	32,85	7,75
40024	LUM	IBL	198,00	33,00	11,00
40027	ELISA	IBL	251,05	49,34	20,68
40029	ELISA	Demeditec	211,56	55,21	35,31
40038	LUM	IBL	191,81	38,20	12,07
40041	ELISA	Demeditec	260,71	62,65	36,18
40049	ELISA	IBL	216,84	54,74	15,55
40055	ELISA	IBL	218,00	51,40	24,70

Sample 1							
Method Supplier	all	LUM IBL	ELISA IBL	ELISA Demeditec	ELISA DRG	ELISA Others	GC/LC-MS Others
Mean	217,4	194,3	241,0	232,9	244,0	221,0	177,4
SD	40,9	12,1	36,7	21,1	69,6		41,5
Mean - 2SD	135,5	170,2	167,7	190,7	104,9		94,4
Mean + 2SD	299,3	218,5	314,4	275,2	383,1		260,3
Min.	117,0	178,0	192,4	211,6	194,8	221,0	117,0
Max.	306,9	215,0	306,9	260,7	293,2	221,0	227,0
Amount	30	6	11	4	2	1	6

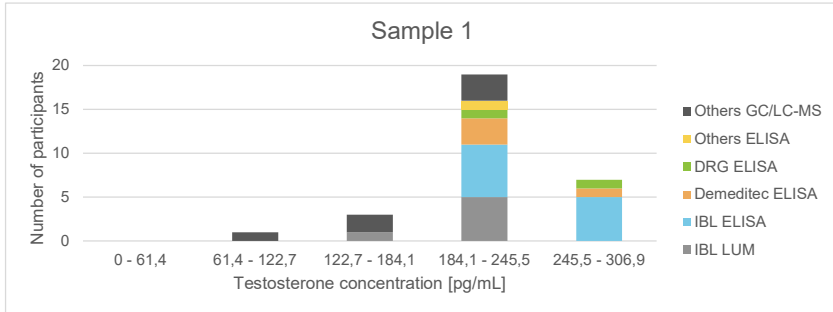

Sample 2							
Method Supplier	all	LUM IBL	ELISA IBL	ELISA Demeditec	ELISA DRG	ELISA Others	GC/LC-MS Others
Mean	45,29	38,29	53,09	61,34	24,86	61,60	33,37
SD	16,08	4,15	15,86	4,12	33,15		5,41
Mean - 2SD	13,12	29,99	21,36	53,10	0,00		22,55
Mean + 2SD	77,45	46,58	84,82	69,58	91,17		44,18
Min.	1,41	33,00	34,18	55,21	1,41	61,60	25,94
Max.	81,18	43,80	81,18	63,90	48,30	61,60	42,70
Amount	31	6	11	4	2	1	7

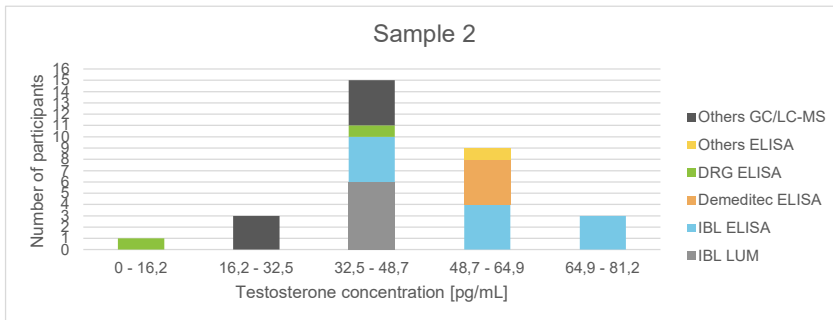

Sample 3							
Method Supplier	all	LUM IBL	ELISA IBL	ELISA Demeditec	ELISA DRG	ELISA Others	GC/LC-MS Others
Mean	18,45	11,29	23,05	37,07	6,86	33,70	7,85
SD	11,80	2,22	9,85	2,51	7,70		2,10
Mean - 2SD	0,00	6,85	3,34	32,05	0,00		3,65
Mean + 2SD	42,05	15,74	42,76	42,10	22,25		12,06
Min.	1,41	7,78	13,94	35,31	1,41	33,70	4,90
Max.	44,81	14,70	44,81	40,80	12,30	33,70	9,90
Amount	31	6	11	4	2	1	7

