



Progesterone			in pg/mL□			
Lab Code	Method	Producer	Sample 1	Sample 2	Sample 3	
3013	ELISA	IBL	63,20	35,70	175,50	
3036	ELISA	DRG	37,41	[1,56]	317,19	
3047	GC/LC-MS	Others	45,87	6,06	158,80	
3050	ELISA	IBL	91,50	40,58	259,11	
3061	ELISA	IBL	78,60	42,10	242,00	
3065	ELISA	IBL	74,69	28,79	198,48	
13004	LUM	IBL	70,30	14,20	184,70	
13009	ELISA	IBL	148,76	63,82	496,37	
13010	ELISA	Salimetrics	85,20	9,40	253,00	
15033	GC/LC-MS	Others	71,00	6,10	253,10	
20002	ELISA	IBL	100,05	32,08	281,75	
20003	ELISA	Others	70,10	21,60	188,30	
30001	LUM	IBL	64,35	17,03	251,74	
40001	ELISA	IBL	95,00	45,00	282,00	
40002	ELISA	IBL	102,10	39,77	259,43	
40003	GC/LC-MS	Others	40,24	[2,70]	135,23	
40014	GC/LC-MS	Others	66,00	[0,00]	199,00	
40015	ELISA	Demeditec	100,20	16,80	270,20	
40016	ELISA	DRG	112,00	21,90	299,70	
40022	GC/LC-MS	Others	47,76	[4,64]	166,21	
40024	ELISA	IBL	89,00	37,00	300,00	
40027	ELISA	IBL	105,05	37,18	242,28	
40049	ELISA	IBL	87,04	35,82	244,76	
40055	ELISA	IBL	66,30	27,00	207,00	
40057	ELISA	IBL	72,00	42,70	248,70	

Sample 1								
Method		LUM	ELISA	ELISA	ELISA	ELISA	ELISA	GC/LC-MS
Supplier	all	IBL	IBL	Salimetrics	Demeditec	DRG	Others	Others
Mean	79,35	67,33	90,25	85,20	100,20	74,71	70,10	54,17
SD	24,54	4,20	22,23			52,74		13,48
Mean - 2SD	30,26	58,92	45,79			0,00		27,21
Mean + 2SD	128,43	75,74	134,72			180,19		81,14
Min.	37,41	64,35	63,20	85,20	100,20	37,41	70,10	40,24
Max.	148,76	70,30	148,76	85,20	100,20	112,00	70,10	71,00
Amount	25	2	13	1	1	2	1	5

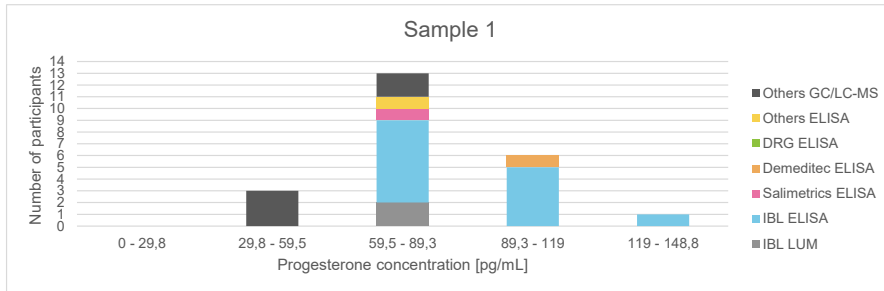

Sample 2								
Method		LUM	ELISA	ELISA	ELISA	ELISA	ELISA	GC/LC-MS
Supplier	all	IBL	IBL	Salimetrics	Demeditec	DRG	Others	Others
Mean	29,55	15,62	39,04	9,40	16,80	21,90	21,60	6,08
SD	14,41	2,00	9,14					0,03
Mean - 2SD	0,74	11,61	20,76					6,02
Mean + 2SD	58,37	19,63	57,32					6,14
Min.	6,06	14,20	27,00	9,40	16,80	21,90	21,60	6,06
Max.	63,82	17,03	63,82	9,40	16,80	21,90	21,60	6,10
Amount	25	2	13	1	1	2	1	5

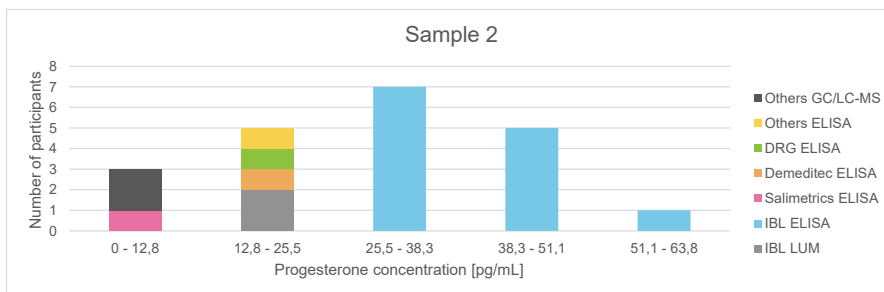

Sample 3								
Method		LUM	ELISA	ELISA	ELISA	ELISA	ELISA	GC/LC-MS
Supplier	all	IBL	IBL	Salimetrics	Demeditec	DRG	Others	Others
Mean	244,58	218,22	264,41	253,00	270,20	308,45	188,30	182,47
SD	69,71	47,40	77,95			12,37		45,60
Mean - 2SD	105,17	123,42	108,52			283,71		91,27
Mean + 2SD	383,99	313,02	420,31			333,18		273,67
Min.	135,23	184,70	175,50	253,00	270,20	299,70	188,30	135,23
Max.	496,37	251,74	496,37	253,00	270,20	317,19	188,30	253,10
Amount	25	2	13	1	1	2	1	5

