



Testosterone

Lab Code	Method	Producer	in pg/mL		
			Sample 1	Sample 2	Sample 3
3013	ELISA	IBL	162,60	25,20	14,10
3019	GC/LC-MS	Others	[862,2]	34,87	8,65
3030	GC/LC-MS	Others	[31,64]	[7,26]	8,10
3036	ELISA	DRG	232,10	35,42	2,65
3047	GC/LC-MS	Others	189,63	46,97	10,09
3050	ELISA	IBL	194,62	27,28	10,88
3061	ELISA	IBL	174,60	39,70	17,50
3065	ELISA	IBL	232,45	45,37	20,61
13004	LUM	IBL	188,80	31,90	7,20
13009	ELISA	IBL	281,14	73,77	34,83
13010	ELISA	DRG	235,20	29,90	4,20
15033	GC/LC-MS	Others	197,70	42,10	9,30
20002	ELISA	IBL	210,39	47,68	12,40
20003	ELISA	Others	202,60	57,80	28,00
30001	LUM	IBL	229,39	37,46	14,70
40001	LUM	IBL	211,00	40,00	17,00
40002	ELISA	IBL	189,52	34,91	11,52
40003	GC/LC-MS	Others	179,56	37,96	7,95
40014	GC/LC-MS	Others	160,00	34,00	6,00
40015	ELISA	Demeditec	222,90	67,20	37,70
40016	ELISA	Demeditec	191,00	54,30	31,20
40022	GC/LC-MS	Others	161,21	34,84	7,64
40024	LUM	IBL	239,00	28,00	12,00
40026	ELISA	Demeditec	177,89	53,21	31,64
40027	ELISA	IBL	214,08	45,96	19,37
40049	ELISA	IBL	191,65	44,65	14,68
40055	ELISA	IBL	223,00	48,10	22,80
40057	ELISA	IBL	175,80	44,00	14,30

Sample 1							
Method		LUM	ELISA	ELISA	ELISA	ELISA	GC/LC-MS
Supplier	all	IBL	IBL	Demeditec	DRG	Others	Others
Mean	202,6	217,0	204,5	197,3	233,7	202,6	177,6
SD	28,3	22,1	33,3	23,1	2,2		16,8
Mean - 2SD	146,1	172,8	137,9	151,0	229,3		144,0
Mean + 2SD	259,1	261,3	271,1	243,6	238,0		211,2
Min.	160,0	188,8	162,6	177,9	232,1	202,6	160,0
Max.	281,1	239,0	281,1	222,9	235,2	202,6	197,7
Amount	28	4	11	3	2	1	7

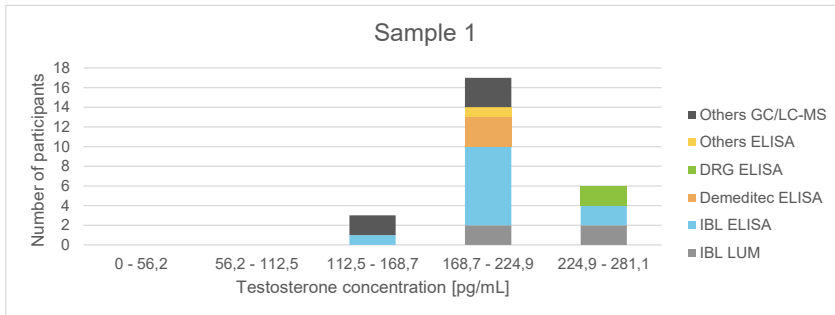

Sample 2							
Method		LUM	ELISA	ELISA	ELISA	ELISA	GC/LC-MS
Supplier	all	IBL	IBL	Demeditec	DRG	Others	Others
Mean	42,32	34,34	43,33	58,24	32,66	57,80	38,46
SD	11,47	5,41	12,85	7,78	3,90		5,14
Mean - 2SD	19,39	23,51	17,63	42,67	24,86		28,19
Mean + 2SD	65,25	45,17	69,03	73,80	40,46		48,73
Min.	25,20	28,00	25,20	53,21	29,90	57,80	34,00
Max.	73,77	40,00	73,77	67,20	35,42	57,80	46,97
Amount	28	4	11	3	2	1	7

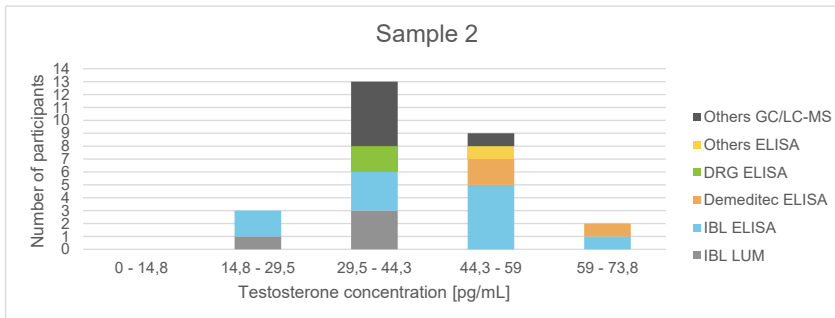

Sample 3							
Method		LUM	ELISA	ELISA	ELISA	ELISA	GC/LC-MS
Supplier	all	IBL	IBL	Demeditec	DRG	Others	Others
Mean	15,61	12,72	17,54	33,51	3,43	28,00	8,25
SD	9,32	4,21	6,91	3,63	1,10		1,30
Mean - 2SD	0,00	4,30	3,73	26,25	1,24		5,64
Mean + 2SD	34,25	21,15	31,35	40,78	5,62		10,85
Min.	2,65	7,20	10,88	31,20	2,65	28,00	6,00
Max.	37,70	17,00	34,83	37,70	4,20	28,00	10,09
Amount	28	4	11	3	2	1	7

